

LABNET INTERNATIONAL, INC.

BioPette™ PLUS

Labnet
Labnet International, Inc.



- Ergonomic design for exceptional comfort
- Dual volume adjustment with volume lock
- Universal tip compatibility
- Fully autoclavable




BioPette™ PLUS

Single Channel Pipette

Labnet's new BioPette Plus pipettes have been engineered to provide the highest level of quality, comfort, accuracy, and precision. Available in single, eight, and twelve-channel pipettes, they provide a range from 0.1µl to 10,000µl.

Each unit is supplied with a calibration certification and a three-year warranty. The BioPette Plus is completely autoclavable and UV resistant. The BioPette Plus chemical-resistant shafts have been designed for universal tip compatibility. The single channel has stainless steel tip ejectors which are easily removed and adjusted to fit any tip.

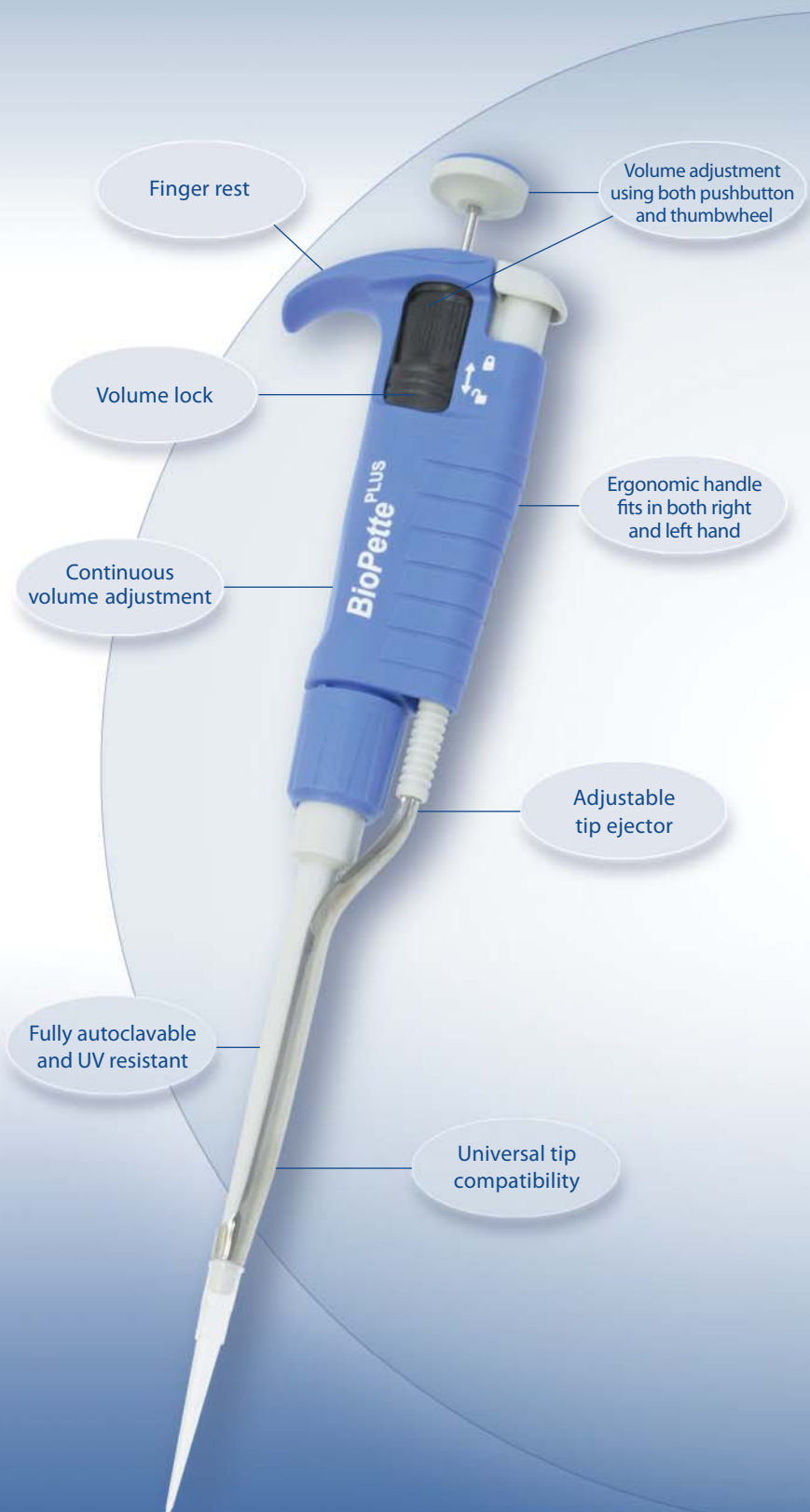

BioPette™ PLUS

Multichannel Pipette

The BioPette Plus finger rest is ergonomically designed, allowing easy single-handed operation and fits comfortably in both right and left hands. The BioPette Plus also features a continuous volume adjustment, and a uniquely designed pushbutton and volume lock which prevents accidental volume change.

While having remarkably low plunger and volume change forces, BioPette Plus pipettes are highly accurate and precise. Each pipette is provided with calibration key for easy in-lab calibration.

BioPette Plus multichannel pipettes feature individually spring-loaded piston assemblies for positive tip sealing. They are provided with a curved ejector bar to produce low ejector forces, and a manifold that rotates 360° for ejecting and pipetting in any direction.



Specifications and Ordering Information

Catalog No.	Volume range	Channels	Accuracy	Precision	Tips
SINGLE CHANNEL PIPETTES*					
P3940- 2	0.1-2µl	1	±12 to ±1.5%	≤6.0 to ≤0.7%	10µl
P3940-10	0.5-10µl	1	±4.0 to ±0.5%	≤4.0 to ≤0.4%	10µl
P3940-20	2-20µl	1	±4.0 to ±0.8%	≤3.0 to ≤0.4%	200µl
P3940-100	10-100µl	1	±1.6 to ±0.8%	≤0.8 to ≤0.3%	200µl
P3940-200	20-200µl	1	±1.2 to ±0.6%	≤0.6 to ≤0.3%	200µl
P3940-1000	100-1000µl	1	±1.0 to ±0.6%	≤0.45 to ≤0.25%	1000µl
P3940-5000	1000-5000µl	1	±0.6 to ±0.5%	≤0.3 to ≤0.2%	5000µl
P3940-10000	1000-10000µl	1	±2.5 to ±0.5%	≤0.6 to ≤0.2%	10000µl
MULTICHANNEL PIPETTES*					
P4508-10	1-10µl	8	±8.0 to ±2.0%	≤6.0 to ≤1.2%	10µl
P4508-50	5-50µl	8	±4.0 to ±1.6%	≤2.5 to ≤0.6%	200µl
P4508-200	20-200µl	8	±3.0 to ±1.0%	≤1.5 to ≤0.6%	200µl
P4508-300	50-300µl	8	±1.6 to ±1.0%	≤1.5 to ≤0.6%	200/300µl
P4512-10	1-10µl	12	±8.0 to ±2.0%	≤6.0 to ≤1.2%	10µl
P4512-50	5-50µl	12	±4.0 to ±1.6%	≤2.5 to ≤0.6%	200µl
P4512-200	20-200µl	12	±3.0 to ±1.0%	≤1.5 to ≤0.6%	200µl
P4512-300	50-300µl	12	±1.6 to ±1.0%	≤1.5 to ≤0.6%	200/300µl
ACCESSORIES					
P3985	Carousel stand for six single channel pipettes				
P3989	Linear rack for six single or multichannel pipettes				
P5423	Multichannel stand holds one 8 or 12 channel pipette				
P3940-SC	Shelf clip holds one single or multichannel pipette				
P3986	Shelf clip holds three single or multichannel pipettes				
P1000-10R	Pipette Tips, 10µl rack of 96				
P1000-10X-R	Pipette Tips, extended length, 10µl rack of 96				
P1000-200R	Pipette Tips, 200µl rack of 96				
P1000-300R	Pipette Tips, 300µl rack of 96				
P1000-1000R	Pipette Tips, 1000µl rack of 100				
P1000-1200R	Pipette Tips, 1200µl rack of 96				
P3998	Pipette Tips, natural, 5000µl, bulk pack 300 tips				
K943	Pipette Tips, clear, 10000µl, bulk pack of 100 tips				
P3940-SK4	BioPette Plus Four Pack includes four pipettes: 0.5 to 10µl , 2.0 to 20µl , 20 to 200µl, 100 to 1000µl Six position carousel stand Three racks of Labnet BioFree™ Pipette Tips				

BioPette Plus Four Pack

