



Conical Competence

The new Eppendorf Conical Tubes 15 mL and 50 mL

The Eppendorf Tubes® – Best Performance for Your Samples

The Eppendorf Conical Tubes 15 mL and 50 mL extend the volume range of Eppendorf's microcentrifuge tubes. They provide life science laboratories with the same outstanding performance for which Eppendorf Tubes are globally recognized as a premium standard.


Safe design and easy handling ensure best sample processing. The Eppendorf Conical Tubes 15 mL and 50 mL are the optimal solution for cell culture and cell biology applications as well as for sample preparation protocols in microbiology and molecular biology laboratories.



Cap with flattened and grooved sides
 > supports a user-friendly ergonomic handling
 > enables a stable upright positioning
 > minimizes contamination risk

Reliable and safe labelling
 due to a large writing area on the tube and due to a flat and light-colored cap

Precise dimensions
 enable maximum compatibility with centrifuge rotors, mixers and shakers

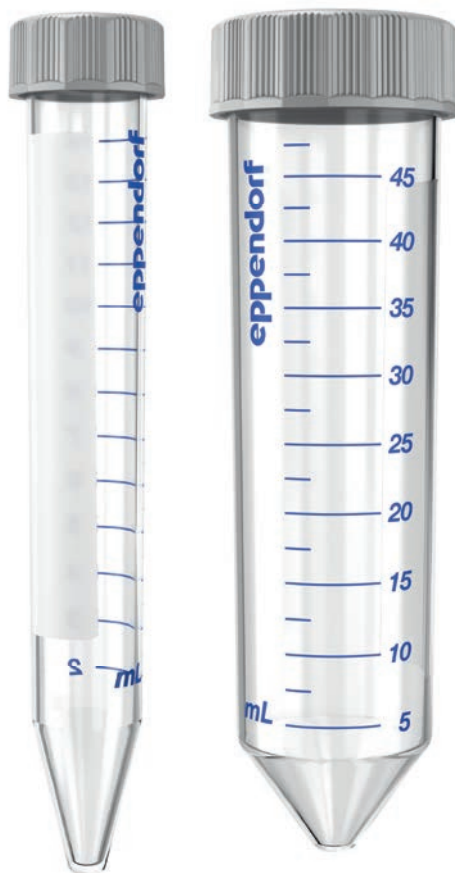

Premium USP class 6 raw material and no use of slip agents, plasticizers, and biocides during manufacturing ensure
 > highest sample integrity
 > high g-Safe® centrifugation stability
 > optimal sample and pellet visibility

Sterile and pyrogen-free
 > DNase-, RNase-, and DNA-free
 > each individual lot is checked and certified for quality and purity



Whether centrifugation, mixing, or other applications — the Eppendorf Conical Tubes 15 mL and 50 mL will perfectly fit into your workflows. They represent the next step in Eppendorf's comprehensive system solutions in liquid, sample and cell handling.

Technical specifications

	Eppendorf Conical Tubes 15 mL	Eppendorf Conical Tubes 50 mL
Material: Tube / Screw cap	Polypropylene (raw material USP class 6) / High-density-Polyethylene (HDPE)	
Dimensions: Height / Diameter	121 mm / 17 mm	116 mm / 30 mm
Operating temperature	-86 °C–100 °C	-86 °C–100 °C
Centrifugation stability	up to 19,500 × g*	up to 19,500 × g*
Certified purity level	sterile, pyrogen-, DNase-, RNase-, and DNA-free	

* Test conditions: Rotor FA-45-6-30 in Eppendorf Centrifuge 5810 R (using a suitable centrifuge adapter for 15 mL tube), sample: Saline solution (density: 1.0 g/mL), controlled centrifugation temp.: 25 °C, centrifugation time: 90 min

Ordering information

Description	International Order no.	North America Order no.
Eppendorf Conical Tubes 15 mL , sterile, pyrogen-, DNase-, RNase-, and DNA-free, 500 Tubes (10 bags of 50 ea.)	0030 122.151	0030122151
Eppendorf Conical Tubes 50 mL , sterile, pyrogen-, DNase-, RNase-, and DNA-free, 500 Tubes (20 bags of 25 ea.)	0030 122.178	0030122178
Eppendorf Tubes®: 0.5 mL – 5.0 mL		
Eppendorf Safe-Lock Tubes, 0.5 mL, Eppendorf Quality™ , colorless, 500 tubes*	0030 121.023	022363611
Eppendorf Safe-Lock Tubes, 0.5 mL, Eppendorf Quality™ , blue, 500 tubes *	0030 121.139	022363646
Protein LoBind Tubes, Protein LoBind, 0.5 mL, PCR clean , 100 tubes (2 bags × 50 tubes)	0030 108.094	022431064
DNA LoBind Tubes, DNA LoBind, 0.5 mL, PCR clean , 250 tubes (5 bags × 50 tubes)	0030 108.035	022431005
Eppendorf Tubes® 3810X, 1.5 mL, Eppendorf Quality™ , colorless, 1000 tubes* ¹	0030 125.150	
Eppendorf Flex-Tube®, 1.5 mL, Eppendorf Quality™ , colorless, 500 tubes* ¹		022364111
Eppendorf Tubes® 3810X, 1.5 mL, Eppendorf Quality™ , blue, 1000 tubes* ¹	0030 125.177	
Eppendorf Flex-Tube®, 1.5 mL, Eppendorf Quality™ , blue, 500 tubes* ¹		022363531
Eppendorf Safe-Lock Tubes, 1.5 mL, Eppendorf Quality™ , colorless, 1000 tubes* ¹	0030 120.086	
Eppendorf Safe-Lock Tubes, 1.5 mL, Eppendorf Quality™ , colorless, 500 tubes* ¹		022363204
Eppendorf Safe-Lock Tubes, 1.5 mL, Eppendorf Quality™ , blue, 1000 tubes* ¹	0030 120.175	
Eppendorf Safe-Lock Tubes, 1.5 mL, Eppendorf Quality™ , blue, 500 tubes* ¹		022363247
Protein LoBind Tubes, Protein LoBind, 1.5 mL, PCR clean , 100 tubes (2 bags × 50 tubes)	0030 108.116	022431081
DNA LoBind Tubes, DNA LoBind, 1.5 mL, PCR clean , 250 tubes (5 bags × 50 tubes)	0030 108.051	022431021
Eppendorf Safe-Lock Tubes, 2.0 mL, Eppendorf Quality™ , colorless, 1000 tubes* ¹	0030 120.094	
Eppendorf Safe-Lock Tubes, 2.0 mL, Eppendorf Quality™ , colorless, 500 tubes* ¹		022363352
Eppendorf Safe-Lock Tubes, 2.0 mL, Eppendorf Quality™ , blue, 1000 tubes* ¹	0030 120.221	
Eppendorf Safe-Lock Tubes, 2.0 mL, Eppendorf Quality™ , blue, 500 tubes* ¹		022363395
Protein LoBind Tubes, Protein LoBind, 2.0 mL, PCR clean , 100 tubes (2 bags × 50 tubes)	0030 108.132	022431102
DNA LoBind Tubes, DNA LoBind, 2.0 mL, PCR clean , 250 tubes (5 bags × 50 tubes)	0030 108.078	022431048
Eppendorf Tubes® 5.0 mL, Eppendorf Quality™ , 200 tubes (2 bags of 100 ea.)	0030 119.401	0030119401
Eppendorf Tubes® 5.0 mL, PCR clean , 200 tubes (2 bags of 100 ea.)	0030 119.460	0030119460
Eppendorf Protein LoBind Tubes 5.0 mL, PCR clean , 100 tubes (2 bags of 50 ea.)	0030 108.302	0030108302
Eppendorf DNA LoBind Tubes 5.0 mL, PCR clean , 200 tubes (4 bags of 50 ea.)* ¹	0030 108.310	0030108310
System accessories		
Rack 5.0 mL (16 spaces) , 2 pcs., white, also suitable for 15 mL conical screw cap tubes	0030 119.495	0030119495
5804/5804R/5810/ 5810 R: rotors and adapters*²		
Rotor FA-45-6-30, for 6 × Eppendorf Conical Tubes 15 and 50 mL	5820 715.006	5820715006
Eppendorf ThermoMixer® C and Eppendorf ThermoStat™ C*²		
Eppendorf SmartBlock™ 15 mL, thermoblock for 8 conical tubes 15 mL	5366 000.021	5366000021
Eppendorf SmartBlock™ 50 mL, thermoblock for 4 conical tubes 50 mL	5365 000.028	5365000028

*¹ Further purity grades, colors and/or package sizes are available at www.eppendorf.com/consumables

*² Further system accessories are available at www.eppendorf.com

Your local distributor: www.eppendorf.com/contact

Eppendorf AG · 22331 Hamburg · Germany
eppendorf@eppendorf.com · www.eppendorf.com
 Eppendorf North America, Inc. · Phone: 800-645-3050
info@eppendorf.com · www.eppendorfna.com
 Eppendorf Canada Ltd. · Phone: 800-263-8715
canada@eppendorf.com · www.eppendorf.ca

www.eppendorf.com