


MTC Bio™

CULTURE & CENTRIFUGE TUBES



- **Round bottom culture tubes**
- **Conical centrifuge tubes**
- **High clarity polymers**
- **Printed graduations & marking spot**
- **Sterile; RNase & DNase free**

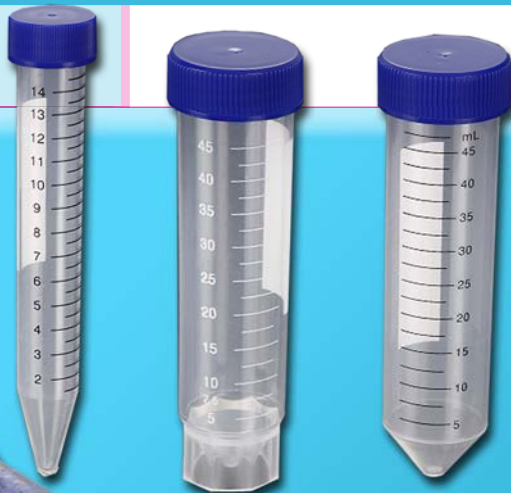

Culture Tubes

These disposable sterile tubes are made of high quality polypropylene (PP) that will withstand temperatures from -190° to +120°C. They are also available in crystal clear polystyrene (PS) with a more limited range of 0° to 70°C. Caps are made of polyethylene (PE) and have the popular two position design (up/ loose for aerobic work and down/ sealed for anaerobic cultures). Polystyrene tubes may be centrifuged at up to 2000xg. Polypropylene tubes withstand higher forces, up to 4000xg.

Centrifuge Tubes

Conical 15ml and 50ml centrifuge tubes are made of premium polypropylene for clarity and centrifugal strength. They are supplied sterile, and although they are autoclavable, their polyethylene screw caps are not. Tubes and caps may be frozen at temperatures as low as -45°C. The maximum centrifugal force for both 15ml and 50ml tubes is 9000xg. Caps have a flat top surface for easy marking and labeling. Cat. No. C2605 features a "skirted" bottom that allows the tube to stand upright without a separate test tube rack.



R1030



R1070

MTC Bio

Div. of MicroSolv Technology Corp.
1 Industrial Way West, Bldg. E, Unit D
Eatontown, NJ 07724 USA
732-578-1777 • www.MTCbio.com

Ordering Information

Culture Tubes (supplied with dual position caps attached)

Cat. No	Capacity	Dimensions	Material	Packaging Unit	Units Per Case	Price Per Case	Qty
T8301	5ml*	12x75mm	PP	25/heat seal bag	40	\$159.00	1000
T8320	5ml*	12x75mm	PP	125/heat seal bag	8	\$99.00	1000
T8300	5ml*	12x75mm	PS	125/heat seal bag	8	\$99.00	1000
T8235	14ml*	17x100mm	PP	25/notched bag	40	\$209.00	1000
T8260	14ml*	17x100mm	PP	1/sterile peel	500	\$129.00	500
T8225	14ml*	17x100mm	PS	25/notched bag	40	\$215.00	1000
T8250	14ml*	17x100mm	PS	1/sterile peel	500	\$129.00	500

* 12x75 tubes are graduated to 4ml. 17x100 tubes are graduated to 12ml

Centrifuge Tubes (supplied with screw caps attached)

Cat. No	Capacity	Dimensions	Material	Packaging Unit	Units Per Case	Price Per Case/500
C2600	15ml	17x118mm	PP	50 in foam rack	10	\$129.00
C2601	15ml	17x118mm	PP	50 in notched bag	10	\$119.00
C2602	50ml	29x115mm	PP	25 in foam rack	20	\$149.00
C2603	50ml	29x115mm	PP	25 in notched bag	20	\$139.00
C2605*	50ml	29x115mm	PP	25 in notched bag	20	\$129.00
R1070	Rack for 30x15ml tubes, pack of 5					\$50.00
R1030	4-way rack for 4 x 50ml, 12 x 15ml, 32 x 1.5ml, 32 x 0.5ml tubes, rainbow pack of 5					\$50.00

* No. 2605 is a self-standing, sterile tube