

STORAGE BOXES

T514 StoreBox™ Storage Boxes

Cover made of polystyrene
Base made of high impact polystyrene

Color your world with a wide variety of economical storage boxes for tubes from 1.2 ml to 10 ml.

These storage boxes are designed to be used at temperatures between -90 °C and +80 °C. Different models are available to accommodate either 25, 42, 81 or 100 sample tubes.

A transparent cover allows the user to see the contents of the box, and is keyed to the base to prevent misalignment. Printed with a series of squares (numbered from 1 to 25, 1 to 42, 1 to 81, or 1 to 100), the surface accepts writing with markers, facilitating inventory control.

A unique color coding system uses colored plastic grids to separate the cover from the base on the 25, 42 and 81-place boxes. Those made to accept 100 tubes (series 2100) have a colored base instead of a grid. Removal of vials facilitated by an innovative vial picker supplied with each storage box (not available with box T514-542). A choice of four popular pastel colors is available. Not autoclavable.




T514-542 StoreBox™ Storage Box

Cover made of polystyrene
Base made of high impact polystyrene

Will hold forty-two 10 ml Sample Tubes. Four colors are available. See page 129 for further details.


Features and benefits of Series 225, 281, 542 & 581 Storage Boxes

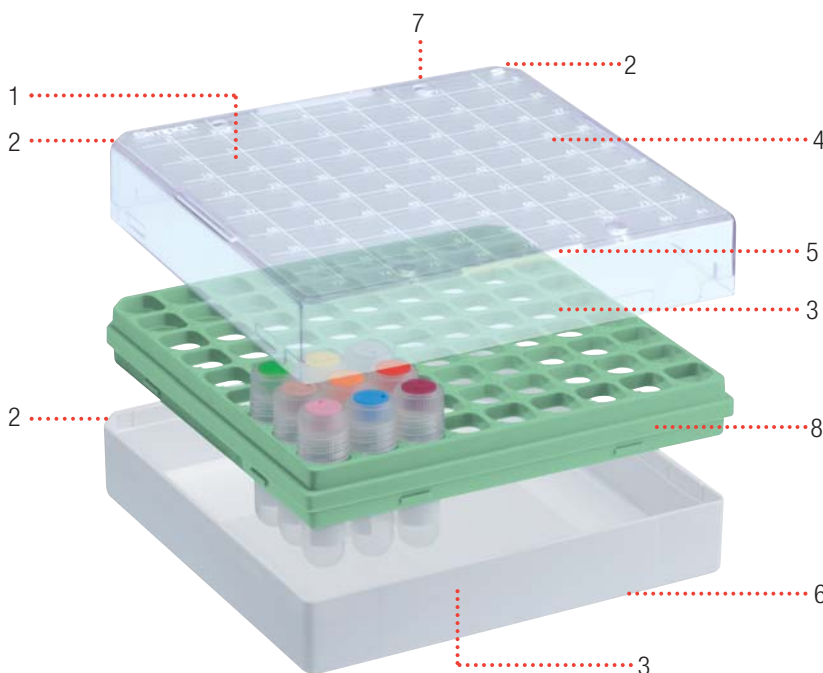

- 1 Cover has numbered squares for easy sample identification
- 2 Two corners of cover and base are cut to prevent misalignment
- 3 Writing surface for identifying base and/or cover
- 4 Vials readily visible through transparent cover
- 5 Air vents minimizing condensation
- 6 Drain holes under base
- 7 Stackable
- 8 Four pastel colors available for better color-coding



A Vial Picker is included with each StoreBox™

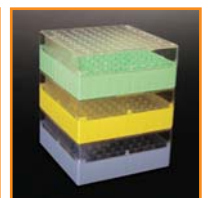

All StoreBox™ Storage Boxes are easily stackable