

# TUBE RACKS

## S500-80 The UniRack™

Made of polypropylene


The UniRack™ offers the laboratory a support far more versatile and easy to use than any other rack available today. It is designed to use minimum counter space while offering maximum flexibility. Made of autoclavable polypropylene, it allows great resistance to various chemicals used in laboratories.

On one side it can hold up to 80 polystyrene or polypropylene 10 and 12 mm tubes, such as 10 x 75 mm or 12 x 75 mm sizes. This rack will accommodate all types of screw cap microtubes from 0.5 to 2 ml made by manufacturers such as Simport, Sarstedt, Nalgene, Bio-Plas, SSI, Sorenson etc... as well as 1 to 5 ml cryogenic vials. Flip the UniRack™ over and you can store up to 60 PCR or microcentrifuge tubes from 0.2 to 0.5 ml.

Units can be firmly anchored laterally to one another, thanks to special anchor pins supplied with each rack. This innovative concept will allow the user to store 80, 160, 240 and even 320 tubes of different shapes, sizes and volumes since the units can be attached to each other either on the 80- or 60- position side facing upward, thus ensuring maximum versatility.

It is supplied with two removable handles allowing for better safety characteristics. The handles of the lower rack make a perfect fit with those of the upper one, ensuring a stable unit which can easily and safely be moved around. An additional protection level is possible by using a very resistant and quite affordable transparent lid allowing a clear view of the contents. More units can be added on top of each other with the lid in place thanks to small pins located on the cover.


There is a frosted area on both sides for bar coding, labeling or writing, enabling the user to identify the contents. It is easy to write on it with a felt-tip pen. These areas are well labeled with arrows to clearly identify which side of the rack the information belongs to. Offered in a wide array of colors allowing the user to classify tubes by their content or by date, work shift, destination or simply by laboratory. Dimensions: 223 x 67 x 27 mm H (9 <sup>3</sup>/<sub>16</sub> x 2 <sup>5</sup>/<sub>8</sub> x 1 <sup>1</sup>/<sub>16</sub> in. H)

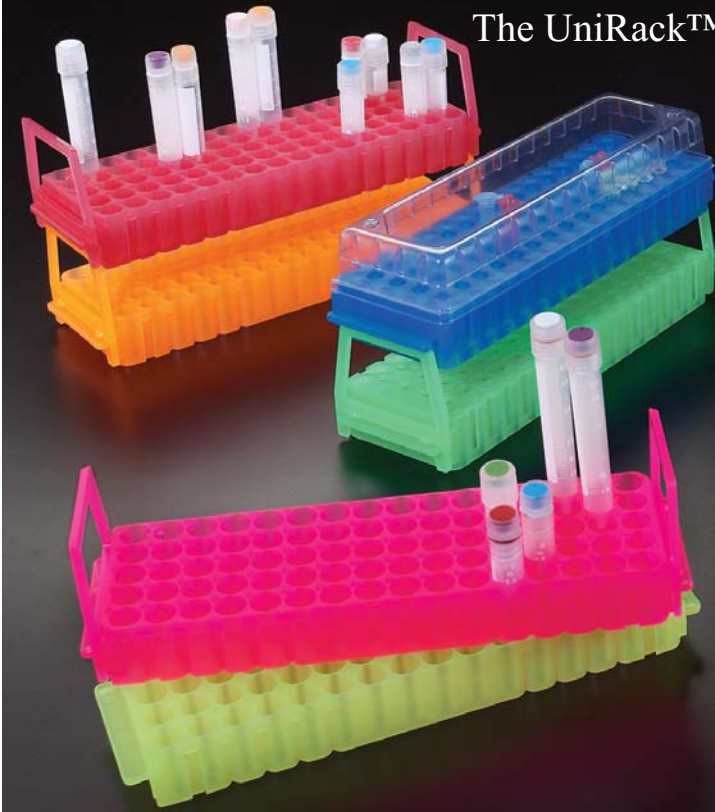


Cat. #	Color	Qty/Cs	Cat. #	Color	Qty/Cs
S500-80B	Blue	10	S500-80R	Red	10
S500-80G	Green	10	S500-80Y	Yellow	10
S500-80O	Orange	10	S500-80AS	Assorted*	10
S500-80P	Pink	10	* Assorted colors : blue, green, orange, pink and yellow		
Lid Cat. #	Color	Qty/Cs			
S501-80	Transparent	10			




The UniRack can also be placed at an angle for easier handling of tubes.

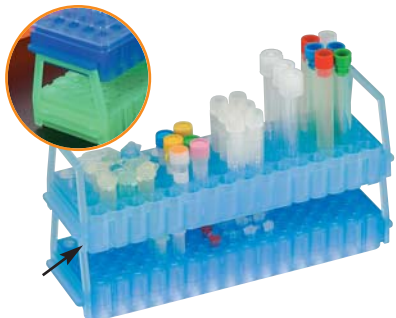

The handles of the lower rack make a perfect fit with those of the upper one, ensuring a stable unit which can easily and safely be moved around.

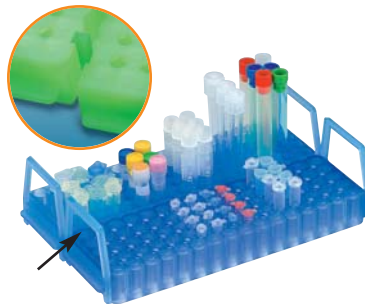

Units can be firmly anchored laterally to one another, thanks to special anchor pins supplied with each rack.

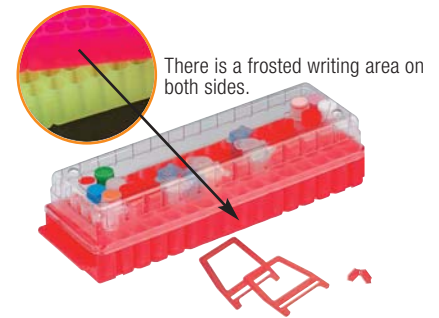

A transparent lid can be placed on the UniRack™, allowing a clear view of the contents.



## S500-25 The UniRack™ Jr.

Made of polypropylene

This smaller model of the UniRack can hold up to 25 polystyrene or polypropylene 10 and 12 mm tubes, such as 10 x 75 mm or 12 x 75 mm sizes. This rack will accommodate all types of screw cap microtubes from 0.5 to 2 ml made by manufacturers such as Simport, Sarstedt, Nalgene, Bio-Plas, SSI, Sorenson etc... as well as 1 to 5 ml cryogenic vials. Flip the UniRack™ over and you can store up to 16 PCR or microcentrifuge tubes from 0.2 to 0.5 ml. Supplied without handles or anchor pins.

Cat. #	Color	Qty/Cs
S500-25B	Blue	10
S500-25R	Red	10
S500-25Y	Yellow	10