


SureSeal™ STERILE SAMPLING BAGS




- **Leak-proof and airtight**
- **Gamma sterilized - no gas residue**
- **Perforated strip - guarantees sterility**
- **Write-on block - for easy labelling**
- **Shielded tabs - no punctures or scratches**

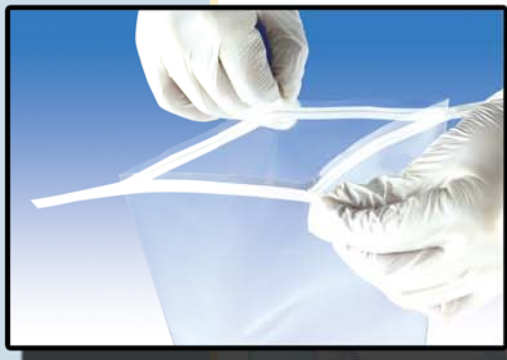

SureSeal™ Bags are gamma sterilized, leaving a virgin interior surface with no toxic residue from gas exposure. Samples can be easily collected by removing the tear strip above the wire closure and pulling on the two external grip tabs to open the mouth of the bag to the desired level. Using the grip tabs ensures sterility and eliminates the risk of contamination from finger contact. Then, introduce the sample into the bag, flip or fold the bag three times, bend the wires and your sample is secure inside the airtight and moisture sealed bag.

SureSeal™ bags are made of virgin polyethylene for excellent clarity and strength. They meet USDA, EPA and FDA standards and are sterile, RNase-free, DNase-free and non-pyrogenic. The PVC tab that surrounds the stainless steel wire is extended beyond the wire ends to eliminate the potential for punctures, scratched hands or ripped gloves due to exposure to uncovered wire ends.

Several sizes are supplied with a white printed area that eliminates tapes, labels, etc. and makes markings more longlasting, clearly visible and less prone to smearing. Although these bags can contain up to 55oz (1600mL), they fold completely flat and require very little storage space.

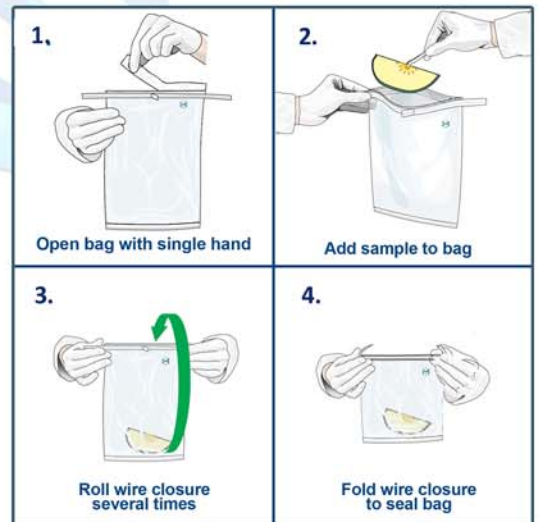

mtc Bio

2600A Main Street Extension
Sayreville, NJ 08872
Ph: 908-279-0280 www.mtcbio.com

Ordering Information

Cat. No.	Description	Capacity	Quantity	Price
B5750	Sterile Sampling Bag, 127mm x 76mm	2oz / 60mL	500	\$58.00
B5760	Sterile Sampling Bag, 178mm x 76mm	4oz / 120mL	500	\$59.00
B5761	Sterile Sampling Bag, 178mm x 76mm, w/ 60x65mm marking spot	4oz / 120mL	500	\$68.00
B5770	Sterile Sampling Bag, 229mm x 114mm	18oz / 530mL	500	\$80.00
B5771	Sterile Sampling Bag, 229mm x 114mm, w/ 60x105mm marking spot	18oz / 530mL	500	\$98.00
B5780	Sterile Sampling Bag, 304mm x 114mm	27oz / 800mL	500	\$115.00
B5790	Sterile Sampling Bag, 381mm x 114mm	36oz / 1060mL	500	\$120.00
B5910	Sterile Sampling Bag, 300mm x 180mm	55oz / 1600mL	250	\$136.00
B5911	Sterile Sampling Bag, 300mm x 180mm, w/ 80x160mm marking spot	55oz / 1600mL	250	\$96.00