

T401 Caps and Stoppers


Made of polyethylene

Plug type push-in caps and 2-position snap caps are made of polyethylene for test tubes with outside diameters of 12 or 17 mm.

Dual position caps offer two possibilities: the closed but unsealed position where samples are maintained aerobic for microbiological procedures; and the fully sealed position where the cap is pushed tight to seal the tube for anaerobic use or for storage, transfer and centrifuge applications. Not autoclavable.



Two-Position Snap Caps

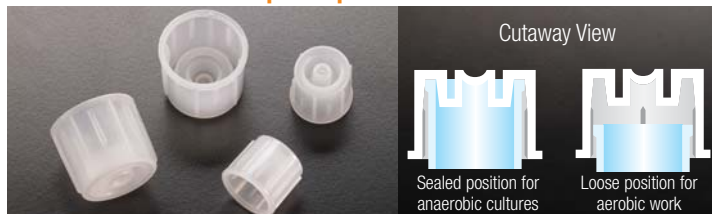

Cat. #	For tubes	Cat. #	For tubes	Color	Qty/Bag	Qty/Cs
T401-3N	12 mm	T401-10N	17 mm	Natural	1000	4000
T401-3B	12 mm	T401-10B	17 mm	Blue	1000	4000
T401-3G	12 mm	T401-10G	17 mm	Green	1000	4000
T401-3R	12 mm	T401-10R	17 mm	Red	1000	4000
T401-3W	12 mm	T401-10W	17 mm	White	1000	4000
T401-3Y	12 mm	T401-10Y	17 mm	Yellow	1000	4000

Cat. #	For tubes made of	For tubes	Color	Qty/Bag
T401-3DSPE	Polystyrene	12 mm	Natural	1000
T401-3DSPP	Polypropylene	12 mm	Natural	1000
T401-10DSPE	Polystyrene	17 mm	Natural	1000
T401-10DSPP	Polypropylene	17 mm	Natural	1000



T402 VACUCAP™ Tube Closures

Made of polyethylene

An economical way to recap blood drawing tubes, disposable glass test tubes and plastic culture tubes. Flexible VACUCAP™ closures protect from aerosols of highly infectious microorganisms. They guard samples against cross-contamination and laboratory work areas against infection and spillage. Precision molded from low-density polyethylene, with a double-flanged seal, VACUCAP™ clamps firmly on the tube. VACUCAP™ holds fast under most rigorous procedures, including centrifugation. Not suitable for autoclaving. Designed for easy-on, easy-off use, due to the exclusive Dual Thumb Tab. Cap removal is simple and quick. VACUCAP™ closure is ideal for recapping most 13 mm and 16mm O.D. evacuated glass blood collection tubes. 13 mm style can also be used on most 12 mm plastic test tubes.

Cat. #	For tubes	Color	Qty/Pk	Qty/Cs
T402-13N	13 mm	Natural	1000	6000
T402-13B	13 mm	Blue	1000	6000
T402-13G	13 mm	Green	1000	6000
T402-13GY	13 mm	Gray	1000	6000
T402-13L	13 mm	Lavender	1000	6000
T402-13R	13 mm	Red	1000	6000
T402-13Y	13 mm	Yellow	1000	6000

Cat. #	For tubes	Color	Qty/Pk	Qty/Cs
T402-16N	16 mm	Natural	1000	6000
T402-16B	16 mm	Blue	1000	6000
T402-16G	16 mm	Green	1000	6000
T402-16GY	16 mm	Gray	1000	6000
T402-16L	16 mm	Lavender	1000	6000
T402-16R	16 mm	Red	1000	6000
T402-16Y	16 mm	Yellow	1000	6000



**A picture is worth a thousand words.
A sample, a thousand pictures...**

You might look at a picture and read the words under it a thousand times, nothing beats having the product in your own hands for evaluation.

Simport is proud to offer you the most comprehensive sample program ever developed in the industry. Just for the asking, you can get free of charge a sample of any Simport product along with a specially designed card describing all the features, benefits and ordering information.

Our Customer Service Specialists are anxiously awaiting your call...
(450) 464-1723